

ZOLL AED 3® Trainer

Training AED and Accessories



Complete ZOLL AED 3 simulation with Real CPR Help real-time CPR feedback

The ZOLL AED 3® Trainer is the only automated external defibrillator (AED) training solution that provides Real CPR Help®, real-time feedback on a trainee's compression quality. The trainer will allow students to learn proper rate and depth of CPR chest compression in a training setting.

The ZOLL AED 3 Trainer is designed to closely mimic the operation of the clinical AED, offering students a more realistic rescue experience. The ZOLL AED 3 Trainer includes:

- An active LCD display with full-colour rescue images that mirror what a student will see on a clinical ZOLL AED 3 device
- Contact sensors in the training pads that detect pad contact with the manikin, eliminating the need for a remote control to support rescue simulations
- An on-device pause button that allows the instructor to pause a training scenario
- Long-lasting consumables, which save time and money
- Easy USB software update/language configuration

2015 AHA Guidelines Ready

Easily modify the configuration and training scenario management on the ZOLL AED 3 Trainer, no remote required. Configurable Options include:


- Training scenario
- Operating mode (Semi-automatic or fully automatic)
- Lay-rescuer prompts
- Check-breathing prompts
- Voice-prompt volume

The ZOLL AED 3 Trainer includes everything an instructor needs:


- ZOLL AED 3 Trainer device
- ZOLL AED 3 Trainer CPR Uni-padz® electrode training harness
- Trainer electrode gels/liner for ZOLL AED 3 Trainer CPR Uni-padz
- Four (4) D cell batteries
- Operator's Guide



Real-time CPR feedback with a large 2-colour compression depth indicator

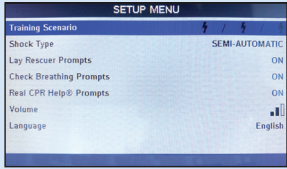

Configure device operation and rescue scenarios directly on the device

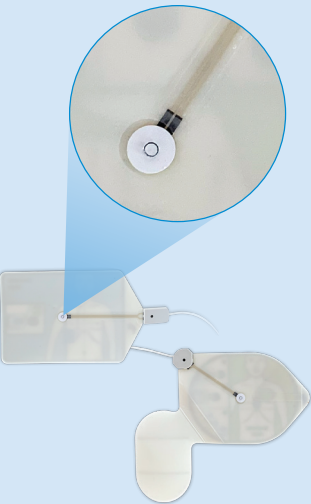

Pad sensors detect contact with the manikin, eliminating the need for a separate remote

ZOLL AED 3 Trainer Specifications

Size: (D x W x H) 5.0 in x 9.3 in x 9.7 in
(12.7 cm x 23.6 cm x 24.7 cm)

Weight: 3 lbs 15oz. (1.8 kg) with batteries

Power: Four (4) D cells (included) for Training AED (Good for more than 16 hours of training)

Part# 8028-000001-01

Rescue Scenarios

Instructors can set up various rescue scenarios directly on the trainer using soft keys.



ZOLL AED 3 Trainer Accessories

Replacements must be ordered separately.



Trainer CPR Uni-padz® electrode replacement adhesive gels
Part# 8028-000014 (case of 5)



Trainer CPR Uni-padz electrode replacement liners Part# 8028-000013



Trainer demo manikin
Part# 8000-0835-01



Trainer USB Software Upgrade Kit
Part # 8028-000011-XX



CPR Uni-padz electrode training harness
Part# 8028-000010

ZOLL Medical Corporation
Worldwide Headquarters
269 Mill Road
Chelmsford, MA 01824
978-421-9655
800-804-4356

For subsidiary addresses and fax numbers, as well as other global locations, please go to www.zoll.com/contacts.

Specifications subject to change without notice.

©2020 ZOLL Medical Corporation. All rights reserved. Real CPR Help, Uni-padz, ZOLL, and ZOLL AED 3, are trademarks or registered trademarks of ZOLL Medical Corporation in the United States and/or other countries. All other trademarks are the property of their respective owners. The ZOLL AED 3 defibrillator is not available for sale in the United States.

Printed in U.S.A. MCN PP 2007 0477

ZOLL